


# Oxygen Flow Splitter

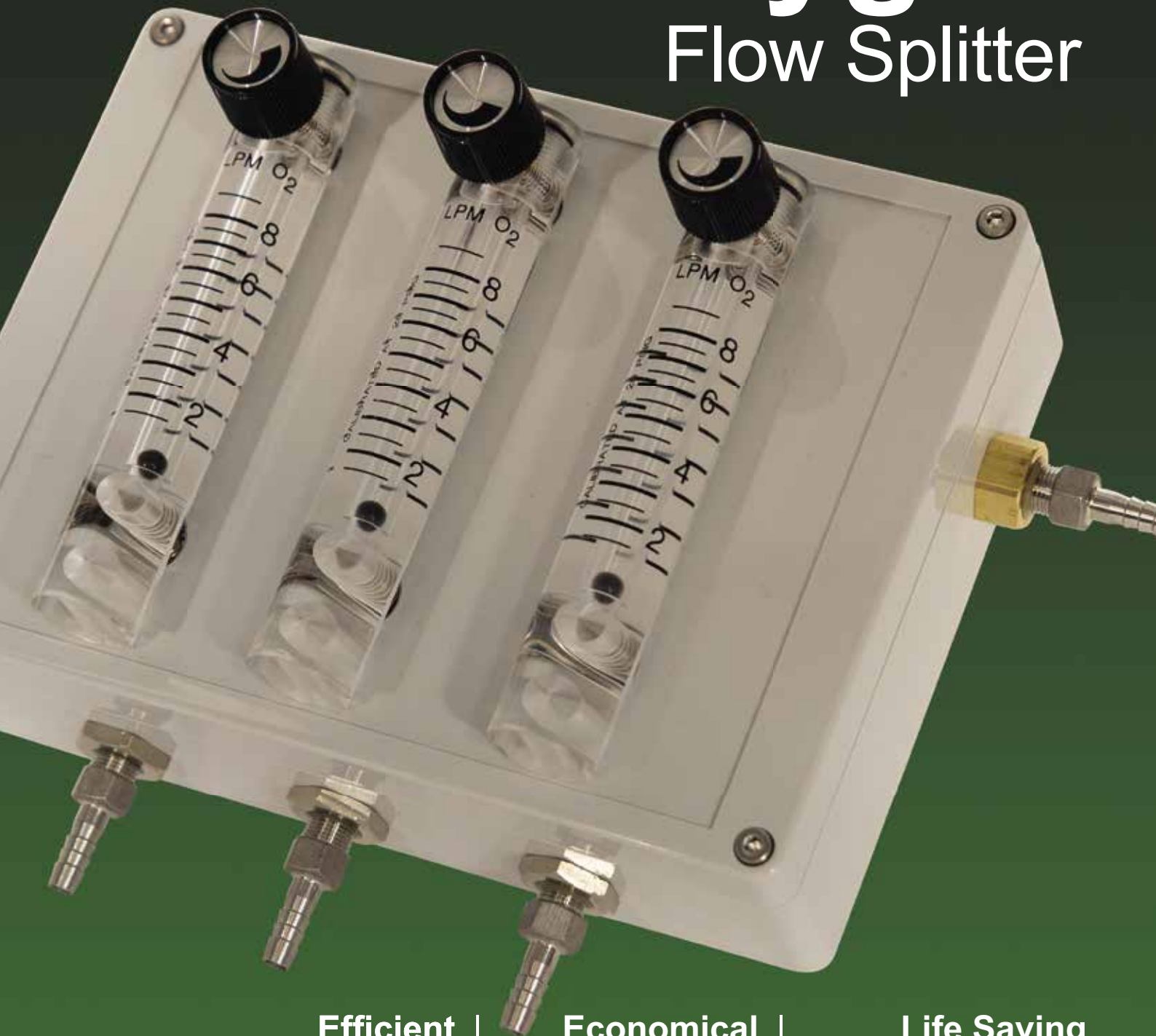

## Efficient

Individually controlled patient oxygen supply from a single source

## Economical

Oxygen titrated to patients at the desired levels

## Life Saving

An essential and effective, oxygen distribution solution

# Oxygen Flow Splitter

Diamedica flow splitters allow the output of a single oxygen source to be split between multiple patients, to provide maximum use of available oxygen without wastage. Once set, split flow rates remain consistent if incoming flow rates increase. Individual adjustable flow meters allow accurate control with a visual indication of the set flow rate.




## Specifications

Patient interface	Compatible with standard oxygen cannula and masks	Flow meter Accuracy	+/- 5% of scale
Flow rate	0-5 / 0-10 lpm per flow (bespoke options available)	Packed Dimensions	16 x 30 x 20 cm, 1.1kg
Inlet Pressure	7-20 psig (48-138 kPa)		

## Key features

- Portable and compact design
- Freestanding or wall-mounted
- Simple user interface
- Multiple patient application
- Simple troubleshooting
- Minimal maintenance
- Minimal cleaning
- No external power required

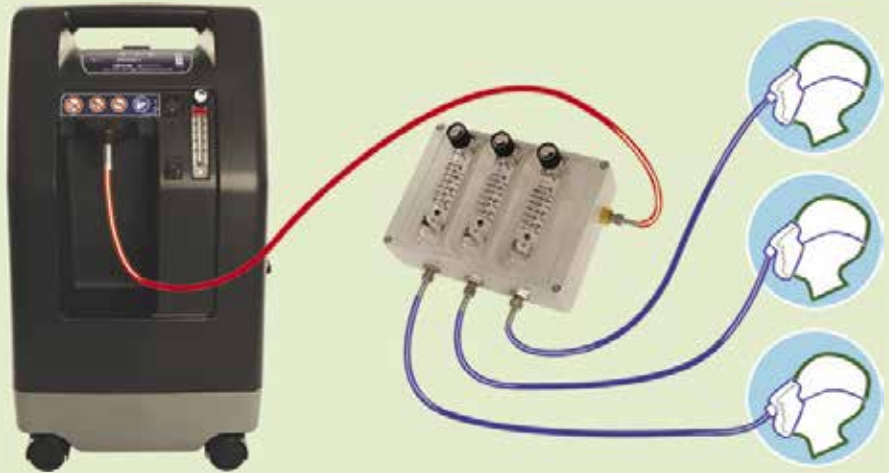

## Training and support

- Full training is available
- Free technical support is available for the lifetime of the machine

 diamedica

**Diamedica (UK) Limited**  
Grange Hill Industrial Estate  
Bratton Fleming  
Barnstaple  
Devon EX31 4UH  
United Kingdom

**Tel:** +44 (0)1598 710066  
**Fax:** +44 (0)1598 710055  
**Email:** [info@diamedica.co.uk](mailto:info@diamedica.co.uk)  
**Web:** [www.diamedica.co.uk](http://www.diamedica.co.uk)

Find us on:    

For free ongoing after sales  
**technical support**



Email: [support@diamedica.co.uk](mailto:support@diamedica.co.uk)



WhatsApp: +44 (0) 7716 503156

© Diamedica (UK) Ltd  
Ref: FLOWSPLIT/2022/01